


## REUSABLE VIDEO LARYNGOSCOPES

by anandic

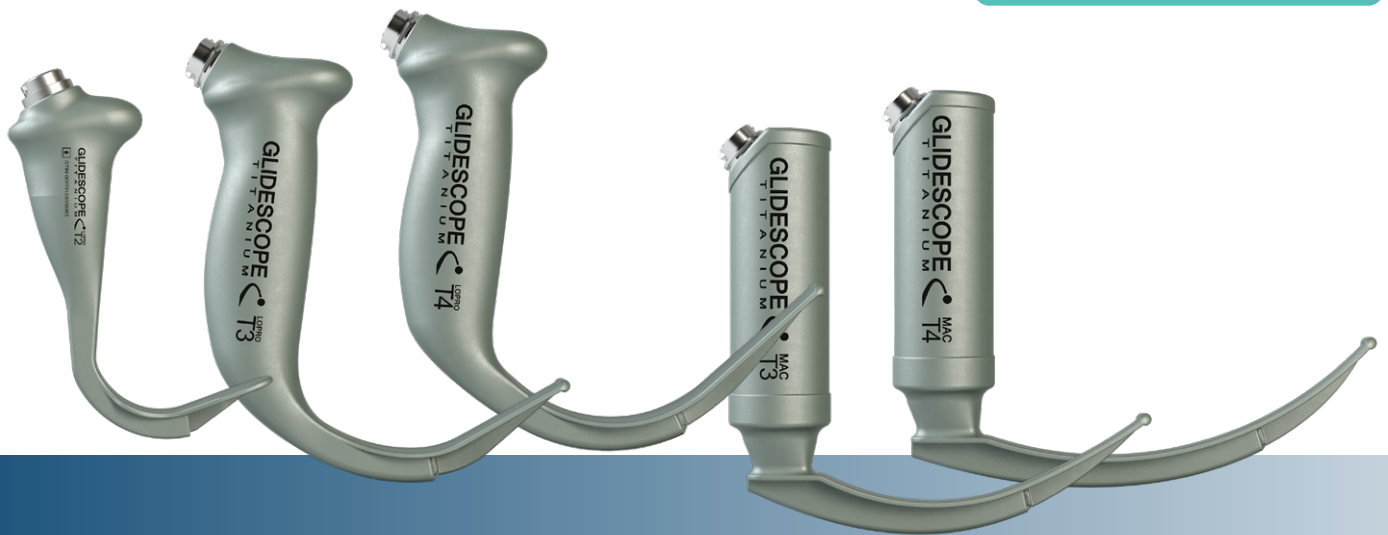

## GLIDESCOPE TITANIUM

The blade design and construction of GlideScope® Titanium video laryngoscopes make them our thinnest, toughest video laryngoscopes ever. These low-profile designs are lightweight and streamlined, offering improved maneuverability and working space for routine and difficult intubations.

GlideScope Titanium video laryngoscopes combine the strength of titanium construction with the confidence of video laryngoscopy. The LoPro T2, T3, and T4 feature our signature GlideScope blade angle for patients from pediatric to large adult. The MAC T3 and T4 provide the familiarity of a Mac blade in a low-profile blade design.

### HIGH-PERFORMANCE FEATURES

- Choice of blade sizes and designs for a wide array of clinical demands
- Low-profile designs for improved maneuverability and working space
- Titanium construction for lightweight durability
- Easy-to-clean blades – IPX8 compliant, with disinfection and sterilization options. No soaking cap required.
- Reveal™ anti-fog technology, with a rapid heating profile to resist lens fogging

### THE TITANIUM SYSTEM

- 5 reusable blade sizes and types
- Compatible with GlideScope® Core™ video monitors
- Mac-style blades for a wide range of clinical needs
- Real-time snapshot, video recording, and video playback features on monitor to help confirm tube placement and facilitate teaching
- External monitor connection capability for education and observation



GlideScope Core 15 with Titanium LoPro T3 and GlideScope Core Video Cable

## APPLICATIONS

GlideScope® Titanium video laryngoscopes provide a clear view of the airway and vocal cords in procedures such as normal or restricted oropharyngeal views/visualization, intubations with blood or other secretions in the airway, airway management in adult and pediatric patients and routine intubations.

## SYSTEM COMPONENTS

- Reusable, titanium blade with digital camera
- 15" or 10" HD touchscreen GlideScope® Core™ monitor
- GlideScope Core workstation with adjustable arm, cable organizer and integrated bins
- GlideScope Core Video Cable for connecting the blade to the monitor
- In-service USB flash drive including the Operations & Maintenance Manual
- 12-volt DC power adapter and power cable
- (Optional) GlideRite® preformed rigid stylets complement angulated video laryngoscopes

## VIDEO LARYNGOSCOPES & ADDITIONAL SYSTEM COMPONENTS

### GlideScope Titanium Video Laryngoscopes



LoPro T2



LoPro T3



LoPro T4



MAC T3



MAC T4

### GlideRite Single-Use Stylets



Small



Medium



Large

### GlideScope Core



GlideScope Core 15 with  
GlideScope Core Premium Workstation

[GlideScope.com](http://GlideScope.com)



**ANANDIC MEDICAL SYSTEMS AG**  
 Stadtweg 24, 8245 Feuerthalen  
[www.anandic.com](http://www.anandic.com)  
[info@anandic.com](mailto:info@anandic.com)  
 Tel. 0848 800 900



[verathon.com](http://verathon.com)

EC REP

Verathon Medical (Europe) B.V.  
 Willem Fenengstraat 13  
 1096 BL Amsterdam  
 The Netherlands  
 Tel: +31 (0) 20 210 30 91  
 Fax: +31 (0) 20 210 30 92